



bed head units

MERYclassic

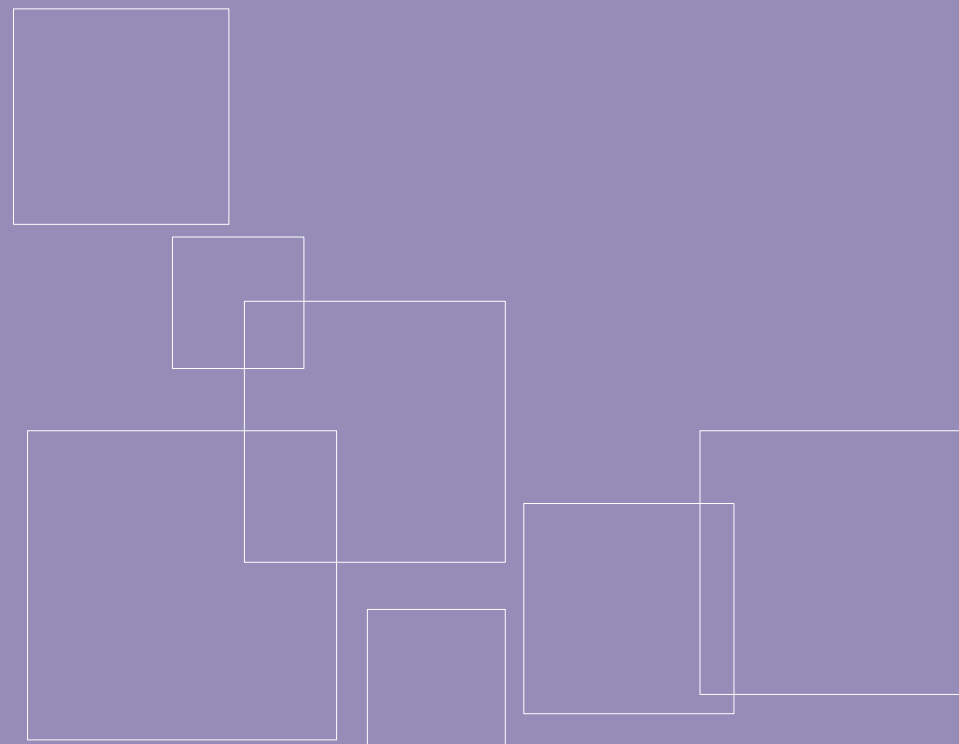

DI: MERY20140325EN

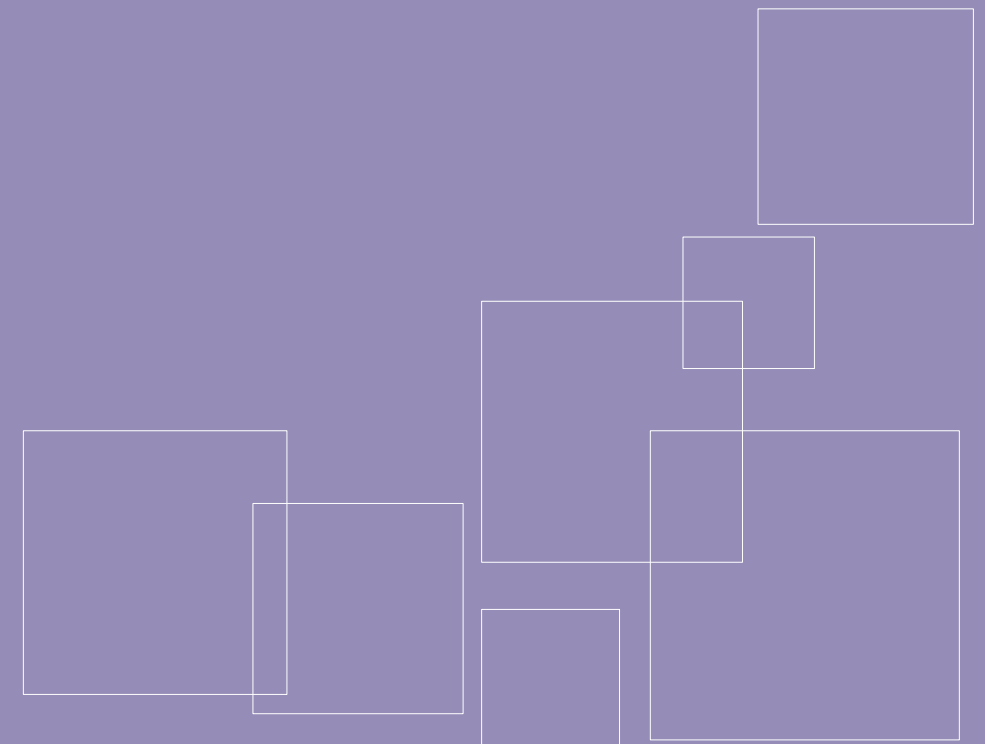



AmbientLight

the classics

MERYclassic is our basic product and the first released in MERY family label. Our focus for the last 15 years has been the constant improvement of its design. In the course of time the design has been modified and improved many times to reach the updated, most modern and most functional form. It has recently reached its perfection and excellence. MERYclassic is the most recognized and appreciated bed head unit by Polish hospitals staff. Its characteristic and outstanding features are highly effective lighting as well as configuration variety and flexibility, granting it the status of leading product in the market. MERYclassic meets the needs of both the clients mostly searching for functionality and the servicers, stressing the importance of the accessibility to media as key issue in maintenance and service with reference to medical gas outlets or electrical sockets.

Spacious internal part and ergonomic management of installation ducts is highly appreciated by fitters. The outstanding features of this bed head unit including its recognition, also expressed by the number of units sold so far, grant MERYclassic the status of leading label in the market.







experience in perfection

MERYclassic is well recognized, ergonomically designed medical device granted IIb class, made of aluminum profiles. Commonly it is selected for standard in-hospital rooms. MERYclassic has unlimited configuration flexibility with reference to the categories of medical gas outlets, IT and electrical sockets including extra equipment like AmbientLight, HandsFree or ControlLight. Ambientlight operates inside standard light cover and through front cover illumination including alternative RGB color selection in real time. HandsFree system may activate indirect and direct light by programmed gestures. ControlLight, on the other hand, controls the light by nurse call system manipulator. MERYclassic provides comfortable room lighting by light reflection and meets all standards applicable. The selection of 5 degrees or 25 degrees angle of upper and lower unit surface where light cover is located, provides accurate light level to the patient. There is milk glass and transparent glass variety available with diversified transparency parameters and frame aluminum reflectors providing most effective light contrasted with minimal energy exhaustion.

The equipment typical for in-hospital rooms:

- standard light T5, 2G11, T8 with single fluorescent light from 13W to 2x80W power ranging or full LED all located in polished asymmetric reflector,
- spotlighting T5, 2G11, T8 from 18W to 36W power ranging or full LED,
- LED night light from 0,6W to 5W power ranging, fluorescent light available,
- gas outlets O2, AIR, VAC based upon detailed customer specification,
- electrical sockets- variety of colors available, modules 45x45 made of or faced with, phase signalization available,
- equipotential bonding pins,
- adjusting to low current installation e.g. RJ45, RJ12 etc.,
- adjusting to any nurse call system,
- 3,4,5, duct profile application including extra equipment fittings e.g. to fit to medical rails etc..



adjustability

Every bed head unit may be multi-positional to integrate medical gas and electrical installation sources in one location. The length of single-positional MERYclassic is 1.6m in compliance with ergonomics and standards applicable. The length of multi-positional unit is a multiple of 1.6m. The length might be adjusted upon request of customer.

Standard version of MERYclassic is equipped with three separate ducts. Two external lighting and electrical canals and an internal gas piping duct, horizontally divided. The clear and easy application of two internal canals, including an open one and an installation between unit and wall makes MERYclassic highly appreciated and recognized by fitters.

