

ultralight

MERYlight

DI: MERY20140325EN



challenge

Hence the real challenge for us so far has been to release a product combining functionality of classic model with modern design. Basing on our technical experience and skills adding a bit of fantasy, we have released an excellently simple and impressively beautiful unit. MERYlight combines the perfect two: functionality and lightness. However its uniqueness is also expressed by its slimness, estheticism and sheer beauty. These features distinguish MERYlight from other products available in the market.

Our engineers and designers employed 9 types of unique aluminum profiles to reach unlimited configuration options. As the result every canal can be fitted with various systemic electrical sockets as well as various gas outlets and slim light covers.

light & lightness


Thanks to combination of uniquely light construction and modern design MERYlight gains rapid recognition in the market. MERYlight is a medical device granted IIb class. This product is made of aluminum and lampshades profiles. MERYlight is characterized by unlimited configuration options with reference to technical and electrical aspects like switching lights by use of nurse call system manipulator or wall switches with optional light-dimming.

The equipment typical for in-hospital rooms:

- standard light T5, 2G11, with single fluorescent light from 21W to 80W power ranging or full LED all located in polished asymmetric reflector,
- spotlight T5, 2G11 standard or from 18W, 2G11 24W, 36W power ranging or full LED,
- LED night light from 0,6W to 5W power ranging, fluorescent light available,
- gas outlets O2, AIR, VAC fitted in main duct,
- electrical sockets- variety of colors and phase signalization available,
- equipotential bonding pins,
- adjusting to low current installation e.g. RJ45, RJ12 etc.,
- multi-module construction equipped with 2, 3 or multi duct profiles including any equipment fittings e.g gas outlets, electrical equipment or any extra equipment to front of unit e.g. to fit to medical rails etc..



Full Light

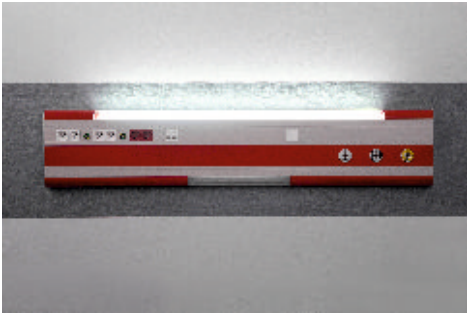

Indirect Light



Direct Light







adjustability

MERYlight is based on aluminum profiles, ensuring the variety of gas outlets, electrical and IT sockets arrangements, fitted directly in ducts to provide clear and simple access, designed solely for this model. The gage of the unit is 65mm. This ultra lightness impacts the cost reduction of reinforcing the wall to which the unit is to be mounted to.

Certainly MERYlight also shares characteristic features of MERY family. It can be multi or single-positional and equipped with extra items like medical appliances carrying systems, shelves, holders etc. *visit MERYcare.

